


Correct packaging

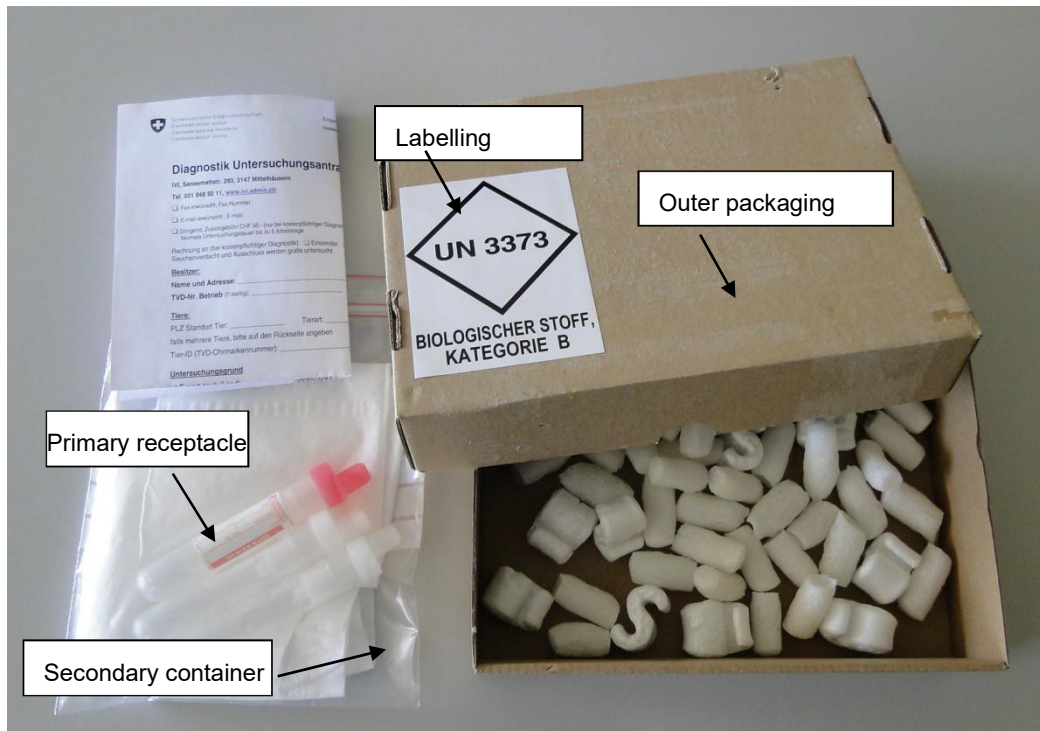

Possible primary receptacles: test tube, leak-proof plastic container

Possible secondary containers: plastic bag, tube, plastic box, sealable canister

Sufficient absorbent material must be added.

Outside packaging: cardboard box, lined envelope

Important: The secondary container or outer packaging must be rigid!

The package must be marked with dangerous goods label UN3373 (for diagnostic material). Label UN3373 can be printed on A4 adhesive paper; printing templates are available at [Diagnostics Mittelhäusern site \(admin.ch\)](http://Diagnostics.Mittelhäusern.site/admin.ch).

Packaging must comply with the following dimensions:

The sides of the UN3373 diamond must be at least 50 mm long. The width of the border around the diamond must be at least 2 mm and all letters and numbers must be at least 6 mm high. The dangerous goods label must be displayed on the package in such a way that the UN 3373 symbol is a diamond (i.e. a square set at 45° as illustrated) (specification according to IATA P650).

(Pre-printed UN3373 dangerous goods labels can also be ordered online from several suppliers).